

BleedSolv™

Stop bleeding fast

Topical Haemostatic Dressings

Product Description

A haemostatic material applied directly to the source of bleeding for use in the control of bleeding

- ❖ REGENERATED CELLULOSE
- ❖ SOLUBLE, WOVEN, TEXTILE
- ❖ pH 6.9

Range of Use

- Exodontia, tissue biopsies, placement of endosseous implants, periodontal surgery, apicoectomy and retrofill, pre-prosthetic surgery, orthognathic surgery and tumour excision.
- Perioperative haemorrhage, intraoperative haemorrhage and oozing or seeping bleeding

Contraindications: The presence of infection. Treat the infection prior to use.

Approvals, Certifications and Standards

Registered with the Therapeutic Goods Administration (TGA)

Rated 91% "better than other haemostatic materials" by Clinicians who reviewed the product for topical use¹

European Certification (CE)

Approved by Food and Drug Administration (FDA) 510k

Advantages

Suitable for use in the presence of coagulation dysfunction, deficiency and disease

No foreign bodies, no allergens, no adverse effects

Soluble – end products: water and polysaccharides

Sterile wrapped and bacteriostatic properties

Cost effective per use

Does not require debridement– soluble, absorbable

pH neutral, no discomfort

Safe to ingest, no taste


Does not require a prescription

Fast, effective control of bleeding

Cost effective per use, less is more. 1-2 x 1.9mm pieces (1 portion) is sufficient for uncomplicated extractions

Product swells when activated. Not for use as a packing material. No need and not cost effective. Packing sockets increases risk of pressure-related complications

Specifications²

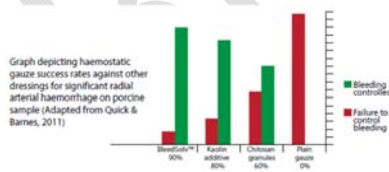

Decreases partial thromboplastin time

Decreases prothrombin release time

Increases platelet and whole blood adhesion rates

Accelerates formation of fibrin cross-linkages

Promotes activation of clotting factor XII

Increases blood viscosity even in the presence of coagulation defects

¹28 Clinicians across 282 uses, 2013

²Against control (plain gauze)

Application	Product Code	Sizes
Exodontia	FADG	19 x 19mm
Procedures, surgery and control of haemorrhage	FAHG: 55,51,11	50 x 50mm, 50 x 100mm, 100 x100mm

